

---

# **BeamScheme Documentation**

***Release 0.52-beta***

**AC Chamberlain**

**Nov 06, 2020**



## CONTENTS:

<b>1</b>	<b>Licence and Disclaimer</b>	<b>1</b>
<b>2</b>	<b>Introduction</b>	<b>3</b>
<b>3</b>	<b>Conventions in this Manual</b>	<b>5</b>
<b>4</b>	<b>System Requirements</b>	<b>7</b>
<b>5</b>	<b>Installation</b>	<b>9</b>
5.1	Uninstalling previous versions . . . . .	9
5.1.1	Under Windows . . . . .	9
5.1.2	Under Linux . . . . .	9
5.2	Installing a new version under Windows . . . . .	9
5.3	Installing a new version under Linux . . . . .	9
<b>6</b>	<b>Running BeamScheme</b>	<b>11</b>
6.1	Running BeamScheme under Windows . . . . .	11
6.2	Running BeamScheme under Linux . . . . .	11
<b>7</b>	<b>Using BeamScheme</b>	<b>13</b>
7.1	The Menu Bar . . . . .	13
7.1.1	The File Menu . . . . .	13
7.1.2	The Export Menu . . . . .	13
7.1.3	The Protocol Menu . . . . .	14
7.1.4	Print PDF . . . . .	16
7.1.5	The Window Menu . . . . .	16
7.1.6	The Help Menu . . . . .	16
7.2	The Profile Panes . . . . .	16
7.2.1	Offset . . . . .	16
7.2.2	Width . . . . .	16
7.2.3	Angle . . . . .	17
7.3	The Image Pane . . . . .	17
7.3.1	Invert . . . . .	17
7.3.2	Normalise to CAX . . . . .	17
7.3.3	Normalise to MAX . . . . .	18
7.3.4	Centre Field . . . . .	18
7.4	The Results Pane . . . . .	18
<b>8</b>	<b>Supported File Formats</b>	<b>19</b>
8.1	DICOM . . . . .	19
8.2	MapCheck . . . . .	19

8.3	PTW	19
8.4	IBA	19
8.5	Image	20
8.6	HIS	20
8.7	XiO	20
8.8	iPlan	20
<b>9</b>	<b>Algorithms</b>	<b>21</b>
9.1	Diagonal/Skew Profiles	21
9.2	Parameters	22
9.2.1	Base Parameters	22
9.2.2	Calculated Parameters	30
9.3	FFF algorithms	38
9.4	In field area	38
<b>10</b>	<b>Protocols</b>	<b>39</b>
10.1	Selecting a protocol	39
10.2	Defined protocols	39
10.2.1	Default	39
10.2.2	AFFSAPS-JORF	40
10.2.3	DIN	40
10.2.4	IEC 60976	41
10.2.5	Elekta	41
10.2.6	Siemens	42
10.2.7	Varian	42
10.2.8	FFF	43
<b>11</b>	<b>Errors and trouble shooting</b>	<b>45</b>
11.1	Unrecognised file or file read error	45
11.2	File error, no data found	45
11.3	This software can not read compressed or 24-bit color files	45
11.4	Unable to load DICOM header segment. Is this really a DICOM compliant file?	45
11.5	This file does not have enough data for the image size	46
11.6	File error, corrupt file	46
11.7	Could not save protocol	46
11.8	No protocol definition files found. Please create a file	46
11.9	Could not read RAW text file	46
11.10	Could not evaluate expression	46
11.11	Unable to fit curve! Initial parameter equal to zero	46
<b>12</b>	<b>Technical Assistance</b>	<b>47</b>
<b>13</b>	<b>Indices and tables</b>	<b>49</b>
	<b>Index</b>	<b>51</b>

## **LICENCE AND DISCLAIMER**

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. The data and formulae on which this program is based are subject to revision and may be changed at any time.
4. You attest that you are competent to understand, evaluate and implement any results generated by the software.
5. Installation, copying or use of this program indicates acceptance of these licence conditions.

Yenzakahle Medical Physics Inc, their subsidiaries and associated companies or practices make no representations or warranties, expressly or implied regarding the program and associated computer files forming the software. No representations or warranties of the merchantability or fitness for any particular purpose are made or that the use of the software will not infringe any patents, copyrights or trademarks. The licensee agrees that neither Yenzakahle Medical Physics Inc nor any of their employees, agents or contractors nor the authors, developers or programmers of the software shall be liable under any claim, charge, demand, whether in contract, tort or otherwise, for any and all loss, cost, charge, claim, demand, fee, expense or damage of any nature or kind arising out of, connected with, resulting from or sustained as a result of this agreement. In no event shall Yenzakahle Medical Physics Inc, nor any of their employees, agents or contractors nor the authors, developers or programmers of the software be liable for special, direct, indirect or consequential damages, losses, costs, charges, demands, fees or expenses of any nature or kind.



## **INTRODUCTION**

Welcome to BeamScheme. This software will assist you in extracting 1D profiles from 2D datasets and calculating various parameters on the profiles. BeamScheme can open various image and 2D array file formats such as SNC MapCheck, PTW 720, IBA Matrix and StarTrack, XiO, DICOM, jpg, etc. Parameters such as field size, field centre, penumbra, flatness and symmetry are calculated. Profiles can be taken at any angle, offset or thickness. Profiles can be exported to a text file for further processing. Results can be printed or exported to PDF. BeamScheme is not intended to replace the commercial software available with 2D arrays, but to complement it.



This manual assumes that you are familiar with using your operating system. If you are not then please consult the appropriate documentation and on-line help. BeamScheme can run on a variety of operating systems and while a thorough knowledge of your specific system is not necessary some basic knowledge such as running files will be an advantage.





## CONVENTIONS IN THIS MANUAL

Various typefaces and punctuation symbols are used to indicate actions that the user must perform. Here is a list of these actions and their meanings.

- Words in a courier font indicate letters that must be typed in on the keyboard. These will usually be data that should be typed into a field.
- Words between single quotes, e.g. 'Start', usually indicate the name of a menu option, field or button.
- Words between angular brackets, e.g. <Enter>, indicate a key that should be pressed. A combination of keys is indicated by a dash '-', for example <Alt-F> means you should press the <Alt> key and the <F> key together. Keys following each other are indicated by a comma ',', for example <Alt-F,S> means that you should press the <Alt> key together with the <F> key, release them and then press the <S> key.
- An icon like this  indicates important information about which you should be aware.
- A hint or short cut is shown like this .

It should be mentioned here that when this manual uses the word 'field' it is referring to a place in the program window where data is entered and not a radiation treatment field unless it is specifically stated.



## **SYSTEM REQUIREMENTS**

BeamScheme can run on a variety of platforms and operating systems. Currently EPS is compiled for:

- Windows 7 or higher
- Fedora 32 running KDE 5 with Qt 5.6

The program needs approximately 5 MB of disk space to install. A 19" monitor or larger is not essential but will make viewing easier. BeamScheme should not be viewed at a resolution of less than 1280x1024. Other linux versions such as GTK based may be available on request.



## INSTALLATION

Installing the Beam Scheme is simple, but if you are concerned that you cannot do it ask your System Administrator or Physicist to do it for you.

### 5.1 Uninstalling previous versions

It is no longer necessary to uninstall previous versions of BeamScheme before installing a new version. Please note that the installation utility will not backup any image or 2D array files in the directory. You will need to do this yourself. If, for any reason, you still want to uninstall the program please follow the instructions below:

#### 5.1.1 Under Windows

In Windows open the 'Control Panel' from the main menu, double click 'Programs and Features', (previously known as 'Add Programs') and click on the 'BeamScheme' entry. Select 'Uninstall' or 'Remove'. When BeamScheme has been removed close 'Programs and Features' and close the 'Control Panel'

#### 5.1.2 Under Linux

Delete the files in the installation directory. If you have any desktop shortcuts or menu links delete those as well.

### 5.2 Installing a new version under Windows

Insert your media containing the installation files. Using Windows Explorer navigate to the drive that you inserted. Select the folder containing the files for your site. Double click on the 'setup' program and follow the on screen instructions. It is best to simply accept the default installation values.

### 5.3 Installing a new version under Linux


Insert the media containing the installation files and mount it. Copy the tar file to the directory where you want to install it. Extract it to this directory. Make sure that the executable has execute permissions. If necessary create a link to the executable and place it in your path.



## **RUNNING BEAMSCHEME**

Contents:

### **6.1 Running BeamScheme under Windows**

Double click the ‘BeamScheme’ icon  on the desktop. Alternatively select ‘Start, All Programs, BeamScheme, BeamScheme’ on the Start Menu.

### **6.2 Running BeamScheme under Linux**

Open a terminal window. Change to the directory where you installed BeamScheme. At the prompt enter:

```
./beamscheme &<enter>
```

Alternatively you can create a link to the application and place it on your desktop. Refer to your system documentation for details.





## USING BEAMSCHEME

Using BeamScheme is simple. There are two profiles, a horizontal or X profile and a vertical or Y profile. The X profile is displayed in the top right quadrant or pane of the screen. The Y profile is displayed in the bottom left quadrant. The result are displayed in the bottom right quadrant. The opened image, fluence or 2D array file is displayed in the top left quadrant. The status bar at the bottom of the screen communicates error messages and warnings. Profiles and results are calculated automatically when a profile is changed or a file is opened.

### 7.1 The Menu Bar

The menu bar allows you to open datasets, export profiles, edit protocols and print results. Some functions are duplicated on the toolbar in *The Image Pane*. Menu options can also be reach using shortcut keys by press <alt> and the underlined character. Available options are:

#### 7.1.1 The File Menu

The file menu options are:

##### File, Open

Opens a dialogue box to select an image, dose map or 2D array file to load. See *Supported File Formats*

This function can also be accessed using the Open File button  on *The Image Pane* toolbar.



If you do not see the file you want make sure you have selected the correct file type.

##### File, Exit

Exits the program



No profiles, data or results are saved.

#### 7.1.2 The Export Menu

The export menu allows you to export the X and Y profiles in text form.

## **Export X Profile**

The X profile can be exported to the clipboard or to a file. The profile is exported as a series of ordered pairs (position, value), with each pair on a new line.

### **Export X profile to file**

Opens a dialogue box to select the location and name of the text file to save the profile to. The profile is saved as a series of ordered pairs (position, value), with each pair on a new line. The file can be imported into a spreadsheet using the CSV format or opened in a text editor.

### **Export X profile to clipboard**

Exports the profile to the clipboard in text format. The profile is saved as a series of ordered pairs (position, value), with each pair on a new line. The data can be pasted into a spreadsheet using the CSV format or opened in a text editor.

## **Export Y Profile**

The Y profile can be exported to the clipboard or to a file. The profile is exported as a series of ordered pairs (position, value), with each pair on a new line.

### **Export Y profile to file**

Opens a dialogue box to select the location and name of the text file to save the profile to. The profile is saved as a series of ordered pairs (position, value), with each pair on a new line. The file can be imported into a spreadsheet using the CSV format or opened in a text editor.

### **Export Y profile to clipboard**

Exports the profile to the clipboard in text format. The profile is saved as a series of ordered pairs (position, value), with each pair on a new line. The data can be pasted into a spreadsheet using the CSV format or opened in a text editor.

## **7.1.3 The Protocol Menu**

The protocol menu allows you to edit protocols and to save them:

### **Editing Protocols**

Selecting “Protocol, Edit” in the menubar allows existing protocols to be edited. The Result window changes to the Edit Protocol window. The parameters are displayed in a grid with the left hand column giving the parameter name and the right hand column giving the parameter equation. Click in the appropriate block to edit. You can also use <ctrl-c> and <ctrl-v> to copy and paste.

To create a parameter section heading leave the equation area blank. Parameter names may be indented by putting spaces in front.

Parameter equations are created by combining base parameters in an arithmetic expression. The expression may not include spaces and should conform to standard pascal.



The multiplication operator is the '\*' not the 'x'.

For example the expression the flatness dose difference paramter

$$100 \cdot \frac{max - min}{max + min}$$

is entered as:

Parameter	Expression
Flatness (%)	100*abs(PMax-PMIn)/(PMax+PMin)

Functions available are:

-  *Save protocol*
-  *Add line*
-  *Delete line*
-  *Exit editing*

## Save protocol

When in editing mode opens the save dialog to save the protocol. It is highly recommended that you do not save over existing protocols. Instead give the changed protocol a unique name. Protocols saved in other directories will not be identified by the program. If you cancel the dialog you will return to the editing mode.




Do not save over a predefined protocol.




Do not save the protocol in any other directory.

## Add line

Inserts a blank line into the parameter list. Click on the line above which you wish to add the new parameter. Click "Add line" . Edit the parameter name and equation.

## Delete line

Deletes a line in the parameter list. Click on the line you wish to delete. Click "Delete line" .

## Exit editing

Aborts editing leaving the parameter set unchanged.

### 7.1.4 Print PDF

From version 0.51 BeamScheme prints the displayed form directly to a pdf. This gives much more flexibility in printing. Selecting the option opens a file dialog for the file name.

If a hard copy is needed the pdf can be printed by any pdf reader.

### 7.1.5 The Window Menu



The Window menu allows you to maximise and/or restore the individual panes as follows:

- Image : Maximise the Image pane
- X Profile : Maximise the X Profile pane
- Y Profile : Maximise the Y profile pane
- Results : Maximise the Results pane
- Restore : Restore all panes to their default configuration

### 7.1.6 The Help Menu

Access this documentation. You may also view an “About” box with a description of BeamScheme, the Licence and a list of credits.

## 7.2 The Profile Panes

The profile pane displays the pixel, dose, fluence or detector values under the corresponding profile on the image pane. The profile pane can be maximised in the application window by clicking the maximise button . Individual profile values can be read by hovering the mouse cursor over the section of profile in question. The profile pane can be returned to its original state by clicking the minimise button .

The profiles are altered by the Offset, Width and Angle spin boxes.

### 7.2.1 Offset

Moves the profile away from the centre of the detector or image. Values can be negative or positive. Values can be incremented or decremented via the up and down arrows in the spin box. Values can be entered manually via keyboard.



Hint: You can change the offset quickly by clicking with the mouse on the image.



### 7.2.2 Width




Widens the profile symmetrically around the profile centre. Values can be increased or decreased by using the up and down arrows respectively in the spin box. Values can be entered by keyboard. Default value is 1. This means the value will always be odd and will increment by 2. Only positive values are allowed. If an even value is entered via the keyboard the profile width will be the next lowest odd value.





### 7.2.3 Angle

Increments or decrements the profile angle. Values can be increased or decreased by using the up and down arrows respectively in the spin box. Values can be entered by keyboard. The zero angle is in the middle of the image left hand size. The X profile starts at angle 0 and the Y profile at angle 90. The X profile can have values between -44 and 45 degrees. The Y profile can have values between 46 and 135 degrees.

## 7.3 The Image Pane

The image pane displays the image, 2D array values, dose or fluence map. The image pane can be maximised in the application window by clicking the maximise button . The image can be windowed using the slider on the right hand side of the image. The image pane can be returned to its original state by clicking the minimise button . You can quickly move the profiles to an appropriate point on the image by clicking on this point with the mouse.

The image pane contains a tool bar with buttons that allow various operations can be carried out on the image. The File Open , Print PDF  and Exit  buttons are shortcuts to the equivalent menu bar functions. Other functions available are:

-  *Invert*
-  *Normalise to CAX*
-  *Normalise to MAX*
-  *Centre Field*

#### 7.3.1 Invert

Inverts the grey scale values using the formula:

$$NZ = max - Z + min$$

Where *max* and *min* are the image maximum and minimum respectively. *Z* is the current pixel value and *NZ* is the new pixel value.



If you don't want an image background level normalise the image first before inverting.

#### 7.3.2 Normalise to CAX

Normalises the image so that the central axis (CAX) value of the image is 100 and the minimum is zero using the formula:

$$NZ = \frac{(Z - min) * 100}{cax - min}$$

Where *cax* and *min* are the image central axis value and minimum respectively. *Z* is the current pixel value and *NZ* is the new pixel value.

### 7.3.3 Normalise to MAX

Normalises the image so that the maximum of the image is 100 and the minimum is zero using the formula:

$$NZ = \frac{(Z - min) * 100}{max - min}$$

Where *max* and *min* are the image maximum and minimum respectively. *Z* is the current pixel value and *NZ* is the new pixel value.

### 7.3.4 Centre Field

Shifts the image by the value of the X and Y field centres. The pixel value at the new index is determined using bi-linear interpolation from the four nearest neighbour pixels from the original image. Image wrap around is not used. Pixels on the leading edge of the image are lost and pixels on the trailing edge are duplicated. This function is intended for fine field adjustments only.



This function changes the original image data.

## 7.4 The Results Pane

Results may be viewed in the Results pane. For a description of the various parameters calculated please see [Parameters](#).

If some parameters are not visible the scroll bar on the right hand side of the pane can be used to scroll down or up to the results.

When the protocol edit mode is selected the Results pane changes to the [Edit Protocol](#) pane.

## SUPPORTED FILE FORMATS

Source	File Extension	Description
<i>DICOM</i>	*.dcm	Dicom CT or RT image files uncompressed
<i>MapCheck</i>	*.txt	MapCheck text files
<i>PTW</i>	*.mcc	PTW 729 array files only at this stage
<i>IBA</i>	*.opg	IBA Matrix and StarTrack opg files
<i>Image</i>	*.bmp	Windows bitmap
	*.tif,*.tiff	Tagged Image format files.
	*.jpg,*.jpeg	Joint Photographic Expert Group files
<i>HIS</i>	*.his	Raw unstructured data
<i>XiO</i>	None	Elekta XiO dose plane export text file
<i>iPlan</i>	None	BrainLab iPlan dose plane text file

As the 2D array files have all been reverse engineered there is absolutely no guarantee of BeamScheme operating correctly on them. If you have a file that is causing BeamScheme to crash or is giving strange results please forward the file along with a description of the problem to [alanphys@yenzakahle.co.za](mailto:alanphys@yenzakahle.co.za).

### 8.1 DICOM

Beamscheme will read most uncompressed DICOM image formats including CT and RI (electronic portal image device). Compressed DICOM files must first be uncompressed using a PACS such as Conquest.

### 8.2 MapCheck

Beamscheme will read SNC MapCheck I and II text files with or without dose calibration.

### 8.3 PTW

Beamscheme will read single plane acquisition PTW 729 array files. 1500 array is not supported at this time.

### 8.4 IBA

Beamscheme will read single plane acquisition IBA StarTrack and Matrixx array files. Multiple frames are summed and averaged

## 8.5 Image

Beamscheme will various 2D image formats such as Windows bitmap (BMP), Tagged Image Format files (TIFF) and Joint Photographic Expert Group (JPEG) files. For colour images only the red channel is imported. It is recommended that only monochromatic images are used for analysis. Colour images can be converted to monochrome using tools such as GIMP or ImageJ.

## 8.6 HIS

Very little information exists on the HIS file format. It seems to have been used by Perkin-Elmer to export X-ray images. Elekta provide it as an image export format. It appears to consist of a 100 byte header (50 16 bit integers) followed by 16 bit integers giving the image data in rows. Integers 9 and 10 appear to give the number of rows and columns respectively.

## 8.7 XiO

The Elekta XiO or Monaco dose plane export is a text file with a 16 line header followed by dose or fluence values. The file is not associated with any particular extension and is identified by the first line value. If dose values were exported in Gy these are converted into cGy on import.



Care should be taken opening XiO files on Windows computers as the CR/LF combination can cause problems.

## 8.8 iPlan

The BrainLab iPlan dose plane export is a text file with a header followed by a separator line followed by one or more dose planes. The dose planes are prefixed by the orientation. The file is identified by the first 8 characters of the first line. If dose values were exported in Gy these are converted into cGy on import.



The dose plane orientation is from left to right and may be mirrored from what is displayed on the iPlan depending on the patient orientation.



## ALGORITHMS

There is an amazing amount of difference in the calculation of supposedly standard parameters between the different vendors. Not only do the definitions of parameters differ but factors such as whether values are interpolated or not or whether field centering is applied have a significant affect on the calculation of the parameters. It is important that the user understand how BeamScheme calculates the various parameters as this will be the first stop in an investigation into any discrepancy between values calculated by BeamScheme and values calculated by vendor supplied software.


BeamScheme makes absolutely no assumptions about the data in contrast to some commercial products which require parameters such as the field size or SSD to be set. It does not perform any beam centering, and no SSD correction is made. No background correction is assumed. The profile or detector centre is defined as half the profile or detector height or width.

The profile is taken from the windowed data. This can have a serious impact on parameters such as field size and penumbra depending on how the window is chosen. The upshot is that defects in the image such as dead pixels or detectors can be reduced by appropriate windowing.

Calculation of the profiles and parameters are discussed in the following sections.

### 9.1 Diagonal/Skew Profiles

Diagonal or skew profiles introduce a number of inaccuracies, which is probably why none of the major vendors offer them. However, if these inaccuracies are understood then skew profiles may yield a lot of useful information. The problem is simple. Consider the following situation.

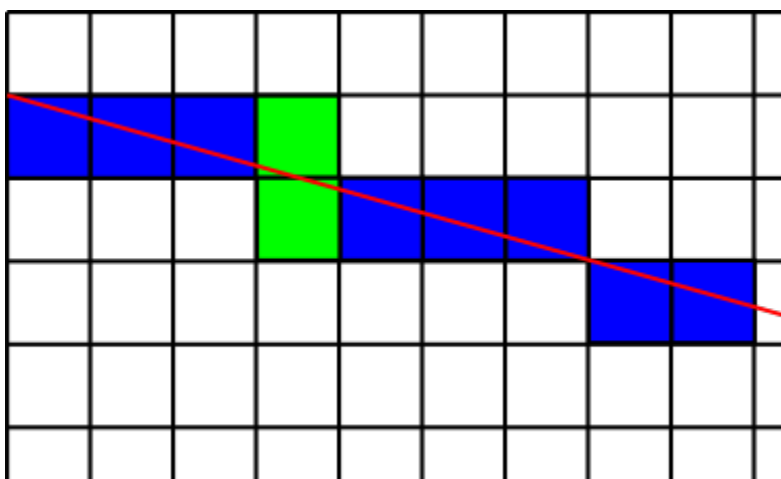

The squares represent pixels or detectors from a 2D dataset. The red line is an ideal skew profile. If the profile completely bisects a square or detector that value is added to the actual profile (blue square). The problem occurs where the ideal profile partially bisects a square or adjacent squares (green squares). There are a number of possible

schemes to solve this such as selecting the nearest pixel or interpolating between adjacent pixels. BeamScheme calculates the length of the profile and then increments an index sequentially along this length. The index is separated into its x and y components which are rounded off to give the pixel indexes along the profile. While the potential does exist for pixels to be skipped or selected twice this is a more robust method and less open to interpretation.

The situation is made worse by the introduction of wide profiles. A wide profile is a profile more than one pixel wide. BeamScheme uses a similar algorithm to that above in that for each index along the profile the normal to the profile is calculated and a second index is sequentially incremented and separated into its x and y components to give the correct pixel. The pixel values are summed. This is done symmetrically on either side of the ideal profile. Thus, the width of the profile will always be odd. This is because an even width will introduce asymmetry in the profile which can lead to further problems.

Skew wide profiles are clipped at the edge of the detector area. This means that the start and finish of wide profiles should be treated with caution and may not be a reflection of reality. The array/image edge will tend to be blurred depending on how skew the profile is.

A similar problem occurs where the profile crosses field edges. If the wide profile is not perpendicular to the field edge this will blur the field edge and increase the penumbra.

## 9.2 Parameters

All parameters are calculated relative to the central axis (CAX) value unless otherwise specified. The CAX is defined as half the profile or detector width or height. Where there is an even number of detectors or pixels and these straddle the CAX the average of the two detectors or pixels is taken.

The parameters, except where otherwise specified, are calculated on the in field area. This area was previously defined as the flattened area, but this definition no longer holds for FFF beams. BeamScheme defines the in field area as 80% of the calculated field size.



Most commercial programs use the nominal field size to calculate the in field area.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are available:

### 9.2.1 Base Parameters

A variety of parameters are automatically calculated on the profiles. While these parameters can be reported directly they are more usually combined in an arithmetic expression. The available base parameters are:

#### Field Statistics

These parameters provide various statistics such as min and max values.



These values are calculated on 80% of the calculated field size.

#### CAX Dose

Dose on the central axis as defined by the array or imaging modality.

Variable name: RCAX



The dose value may be affected by the normalisation or windowing.

### Maximum Dose

Maximum dose of the profile over the in field area.

Variable name: PMax



The dose value may be affected by the normalisation or windowing.

### Minimum Dose

Minimum dose of the profile over the in field area.

Variable name: PMin



The dose value may be affected by the normalisation or windowing.

### Maximum position

Position of the maximum dose point in the profile over the in field area.

Variable name: MPos

### Field Sum

Sum of the dose values of the profile over the in field area.

Variable name: PSum



The dose value may be affected by the normalisation or windowing.

### Field Sum of Squares

Sum of the squared dose values of the profile over the in field area. Can be used to calculate the standard deviation.

Variable name: PSSqr



The dose value may be affected by the normalisation or windowing.

### Number of Profile Points

Number of measurement points in the profile evaluated over the in field area.

Variable name: N



This is not the total number of points in the profile.

## **Interpolated Parameters**

These parameters are the distances from the central axis as defined by the array or imaging modality for a certain dose level and form the base parameters of the traditional definitions of field size and centre. These dose levels usually fall between two measured points. In order to improve accuracy, particularly for 2D arrays where the measurements are relatively sparse, the distance of the dose level is interpolated by linear interpolation from the adjacent data points.

### **Left 90% dose distance**

Distance from the central axis as defined by the array or imaging modality to the point on the left side of the displayed profile where the dose profile drops below 90% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

### **Right 90% dose distance**

Distance from the central axis as defined by the array or imaging modality to the point on the right side of the displayed profile where the dose profile drops below 90% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: R90

### **Left 80% dose distance**

Distance from the central axis as defined by the array or imaging modality to the point on the left side of the displayed profile where the dose profile drops below 80% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: L80

### **Right 80% dose distance**

Distance from the central axis as defined by the array or imaging modality to the point on the right side of the displayed profile where the dose profile drops below 80% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: R80

### **Left edge interpolated**

Distance from the central axis as defined by the array or imaging modality to the point on the left side of the displayed profile where the dose profile drops below 50% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: LEdge

### **Right edge interpolated**

Distance from the central axis as defined by the array or imaging modality to the first point on the right side of the displayed profile where the dose profile drops below 50% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: REdge

### Left 20% dose distance

Distance from the central axis as defined by the array or imaging modality to the point on the left side of the displayed profile where the dose profile drops below 20% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: L20

### Right 20% dose distance

Distance from the central axis as defined by the array or imaging modality to the point on the right side of the displayed profile where the dose profile drops below 20% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: R20

### Left 10% dose distance

Distance from the central axis as defined by the array or imaging modality to the point on the left side of the displayed profile where the dose profile drops below 10% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: L10

### Right 10% dose distance

Distance from the central axis as defined by the array or imaging modality to the point on the right side of the displayed profile where the dose profile drops below 10% of the central axis dose. If this point lies between two measured values the distance is linearly interpolated.

Variable name: R10

## Inflection Point Parameters

For FFF beams the standard interpolated or dose level parameters do not give good results. The field edge of FFF beams is defined at the inflection point of the penumbra. For low resolution measuring devices like 2D arrays this can be quite inaccurate. Therefore a sigmoid model is fitted to the penumbra using a Hill function:

$$f(x) = A + \frac{B - A}{1 + \left(\frac{C}{x}\right)^D}$$

where:

- A: sigmoid low level
- B: sigmoid high level
- C: approximate inflection point
- D: slope of the sigmoid

The inflection point is determined from:

$$x = C \cdot \left( \frac{D-1}{D+1} \right)^{\frac{1}{D}}$$

Once the regression parameters have been determined the inverse Hill function can be used to determine other parameters around the inflection point:

$$x = C \cdot \left( \frac{f(x) - A}{B - f(x)} \right)^{\frac{1}{D}}$$

The profile values to fit are taken from the out field area, that is from 80% of the calculated field size to the end of the profile for fields larger than and including 7x7 cm and are taken from the origin to the end of the profile for fields smaller than 7x7 cm.



If the penumbra is not well formed the non-linear regression will fail and the results returned will be 0.

### Left inflection point

This is the inflection point on the left side of the profile as determined by non-linear regression on a Hill function:

Negative axis increments are converted to positive. The axis sign is then applied to the results.

Variable name: LInf



If the penumbra straddles the origin the non-linear regression will fail.

### Right inflection point

This is the inflection point on the right side of the profile as determined by non-linear regression on a Hill function:

Negative axis increments are converted to positive. The axis sign is then applied to the results.

Variable name: RInf



If the penumbra straddles the origin the non-linear regression will fail.

### 80% left inflection point distance

The traditional parameters for defining the penumbra do not work for FFF beams. The closest comparable value for FFF beams to the 80% dose level distance for WFF beams is to take the distance of 1.6 x the dose value of the inflection point. For WFF beams this corresponds to the 80% dose level (50% x 1.6 = 80%).

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 1.6 x the dose at the inflection point on the left side of the profile.

Variable name: IL80

### 80% right inflection point distance

The traditional parameters for defining the penumbra do not work for FFF beams. The closest comparable value for FFF beams to the 80% dose level distance for WFF beams is to take the distance of 1.6 x the dose value of the inflection point. For WFF beams this corresponds to the 80% dose level ( $50\% \times 1.6 = 80\%$ ).

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 1.6 x the dose at the inflection point on the right side of the profile.

Variable name: IR80

### 50% left inflection point distance

The inflection point is not necessarily at the 50% dose level for WFF beams. To provide a comparative value for WFF beams the distance to the 50% dose level of the Hill function is provided.

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 0.5 x the difference between the sigmoid high level and the sigmoid low level on the left side of the profile.

Variable name: IL50

### 50% right inflection point distance

The inflection point is not necessarily at the 50% dose level for WFF beams. To provide a comparative value for WFF beams the distance to the 50% dose level of the Hill function is provided.

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 0.5 x the difference between the sigmoid high level and the sigmoid low level on the right side of the profile.

Variable name: IR50

### 20% left inflection point distance

The traditional parameters for defining the penumbra do not work for FFF beams. The closest comparable value for FFF beams to the 20% dose level distance for WFF beams is to take the distance of 0.4 x the dose value of the inflection point. For WFF beams this corresponds to the 20% dose level ( $50\% \times 0.4 = 20\%$ ).

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 0.4 x the dose at the inflection point on the left side of the profile.

Variable name: IL20

### 20% right inflection point distance

The traditional parameters for defining the penumbra do not work for FFF beams. The closest comparable value for FFF beams to the 20% dose level distance for WFF beams is to take the distance of 0.4 x the dose value of the inflection point. For WFF beams this corresponds to the 20% dose level ( $50\% \times 0.4 = 20\%$ ).

This is the distance from the origin as determined by the array or imaging modality to the dose level corresponding to 0.4 x the dose at the inflection point on the right side of the profile.

Variable name: IR20

### Inflection point left slope

The slope of the Hill function at the inflection point provides a good indication of the width of the penumbra. This is the slope of the penumbra on the left side of the profile.

Negative axis increments are converted to positive. The axis sign is then applied to the results.

Variable name: LSlope



If the penumbra straddles the origin the non-linear regression will fail.

### Inflection point right slope

The slope of the Hill function at the inflection point provides a good indication of the width of the penumbra. This is the slope of the penumbra on the right side of the profile.

Negative axis increments are converted to positive. The axis sign is then applied to the results.

Variable name: RSlope



If the penumbra straddles the origin the non-linear regression will fail.

### Dose Point Parameters

The traditional flatness parameters have no meaning for FFF beams. Therefore the dose at various points along the profile are used to describe the shape of the profile. The dose values at 20, 50, 60 and 80% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point can be reported.

#### Dose point 20% left

The dose value at 20% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the left side of the profile.

Variable name: LD20

#### Dose point 20% right

The dose value at 20% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the right side of the profile.

Variable name: RD20

#### Dose point 50% left

The dose value at 50% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the left side of the profile.

Variable name: LD50



### Dose point 50% right

The dose value at 50% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the right side of the profile.

Variable name: RD50

### Dose point 60% left

The dose value at 60% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the left side of the profile.

Variable name: LD60

### Dose point 60% right

The dose value at 60% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the right side of the profile.

Variable name: RD60

### Dose point 80% left

The dose value at 80% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the left side of the profile.

Variable name: LD80

### Dose point 80% right

The dose value at 80% of the distance from the origin as defined by the array or imaging modality to the edge of the field defined by the inflection point on the right side of the profile.

Variable name: RD80

## Symmetry Base Parameters

The symmetry parameters apply equally to FFF and WFF fields. These are the base parameters from which the symmetry parameters are calculated.

### Area left

The integration or summation of the left side of the profile from the origin as defined by the array or imaging modality to the edge of the field as defined by the 50% dose value for WFF beams or the inflection point in the case of FFF beams.

Variable name: ALeft



Small sharply peaked fields may be incorrectly identified as FFF beams.

## Area right

The integration or summation of the right side of the profile from the origin as defined by the array or imaging modality to the edge of the field as defined by the 50% dose value for WFF beams or the inflection point in the case of FFF beams.

Variable name: ARight



Small sharply peaked fields may be incorrectly identified as FFF beams.

## Maximum dose ratio base parameter

The maximum dose ratio symmetry (also called the point difference quotient) is the maximum ratio between the left and right profile values at the same distance from the profile or detector centre:


$$\text{maximum} \left[ \frac{P(dL)}{P(dR)}, \frac{P(dR)}{P(dL)} \right]$$

for  $dR = -dL$  from the origin as defined by the array or imaging modality to the end of the in field area.

Variable name: RDiff



The ratio symmetry is affected by the field centre. If the field is slightly offset you can use the “Centre field”

tool  to correct any offset.

## Point difference base parameter

The point difference (also called maximum variation) is the maximum absolute difference between the left and right profile values at the same distance from the profile or detector centre:

$$P(dL) - P(dR)$$

for  $dR = -dL$  from the origin as defined by the array or imaging modality to the edge of the in field area.

Variable name: ADiff



The point difference symmetry may be affected by the field centre. If the field is slightly offset you can use the

“Centre field” tool  to correct any offset.

## 9.2.2 Calculated Parameters

All parameters are calculated relative to the central axis (CAX) value unless otherwise specified. The CAX is defined as half the profile or detector width or height. Where there is an even number of detectors or pixels and these straddle the CAX the average of the two detectors or pixels is taken.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters have been defined in the existing protocols:

## Field Parameters

The field parameters describe the geometry of the field. The edge of the field can be described in a number of ways. The most common are the Full Width Half Maximum (FWHM) or field size at 50% of the CAX value, and the inflection point.

## Interpolated Field Parameters

All parameters are calculated relative to the central axis (CAX) value unless otherwise specified. The CAX is defined as half the profile or detector width or height. Where there is an even number of detectors or pixels and these straddle the CAX the average of the two detectors or pixels is taken.

## Field Edge

The field edges are defined at the first point from the centre of the profile where the profile reaches the 50% value, i.e.:

$$\frac{CAX}{2}$$

Where *CAX* is the value of the profile on the central axis (CAX).

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.

The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. If the 50% value lies between measured points the field edge is interpolated. The value reported is the linear distance from the profile centre to the field edge. The left and right field edge of the profile is given.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to a reduced field edge value being reported.

## Field Centre

The field centre (*FC*) is defined as the distance in cm from the centre of the detector to the centre of the field and is calculated as the half the difference between the distance from the profile centre to the right and left *Field Edge*:

$$FC = \frac{d(RE) + d(LE)}{2}$$

Where *d(RE)* and *d(LE)* are the linear distances from the profile centre to the right and left field edges respectively.

## Field Size

The field size (*FS*) is defined as the distance from the left *Field Edge* to the right *Field Edge*:

$$FS = d(RE) - d(LE)$$

Where *d(RE)* and *d(LE)* the linear distances from the profile centre to the right and left field edges respectively.



Noisy data can lead to an inaccurate field size being reported.



The right and left field edges are affected by the set window levels.

## Inflection Point Field Parameters

All parameters are calculated relative to the inflection point of the penumbra unless otherwise specified. The inflection point is determined by fitting a Hill function to the penumbra data.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are available:

### Inflection Point Field Edge

The field edges as given by the inflection point are defined at the first point from the centre of the profile where the gradient of the profile reaches its maximum value. This is the same as the inflection point of a Hill function fitted to the penumbra.:

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.

Refer to the base *Inflection Point Parameters* for more detail on the fitting of the Hill function.



The profile must be centered on the origin as defined by the array or imaging modality for best results.

### Inflection Point Field Centre

The field centre (*FC*) is defined as the distance in cm from the centre of the detector to the centre of the field and is calculated as the half the difference between the distance from the profile centre to the right and left *Inflection Point Field Edge*:

$$FC = \frac{d(RE) + d(LE)}{2}$$

Where  $d(RE)$  and  $d(LE)$  are the linear distances from the profile centre to the right and left field edges respectively.



The profile must be centered on the origin as defined by the array or imaging modality for best results.

### Inflection Point Field Size

The field size (*FS*) is defined as the distance from the left *Inflection Point Field Edge* to the right *Inflection Point Field Edge*:

$$FS = d(RE) - d(LE)$$

Where  $d(RE)$  and  $d(LE)$  the linear distances from the profile centre to the right and left field edges respectively.



The profile must be centered on the origin as defined by the array or imaging modality for best results.

## Penumbra Parameters

Penumbra parameters have traditionally been defined according to the dose level. Where the dose level falls between measurement points the distance to the dose level is interpolated by linear regression. This approach does not work for FFF beams and instead the penumbra is defined relative to the inflection point of the penumbra.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are defined in the existing protocols:

## Interpolated Penumbra Parameters

All parameters are calculated relative to the central axis (CAX) value unless otherwise specified. The CAX is defined as half the profile or detector width or height. Where there is an even number of detectors or pixels and these straddle the CAX the average of the two detectors or pixels is taken.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are available:

### Penumbra 90%-10%

The 90%-10% penumbra is defined as the distance from the first point on the profile at 90% of the CAX value to the first point that reaches 10% of the CAX value:

$$Penumbra(90 - 10) = abs(d(P_{10\%}) - d(P_{90\%}))$$

$$P_{10\%} = CAX * 0.1$$

$$P_{90\%} = CAX * 0.9$$

Where  $d(P_{10\%})$  and  $d(P_{90\%})$  are the linear distances to the 10% and 90% profile values respectively as defined above. The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. The left and right field penumbra of the profile is given.

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to an inaccurate penumbra value being reported.

### Penumbra 80%-20%

The 80%-20% penumbra is defined as the distance from the first point on the profile at 80% of the CAX value to the first point that reaches 20% of the CAX value:

$$Penumbra(80 - 20) = abs(d(P_{20\%}) - d(P_{80\%}))$$

$$P_{20\%} = CAX * 0.2$$

$$P_{80\%} = CAX * 0.8$$

Where  $d(P_{20\%})$  and  $d(P_{80\%})$  are the linear distances to the 20% and 80% profile values respectively as defined above. The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. The left and right field penumbra of the profile is given.

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to an inaccurate penumbra value being reported.

## Penumbra 90%-50%

The 90%-50% penumbra is defined as the distance from the first point on the profile at 90% of the CAX value to the first point that reaches 50% of the CAX value or field edge:

$$Penumbra(90 - 50) = abs(d(P_{50\%}) - d(P_{90\%}))$$

$$P_{50\%} = CAX * 0.5$$

$$P_{90\%} = CAX * 0.9$$

Where  $d(P_{50\%})$  and  $d(P_{90\%})$  are the linear distances to the 50% and 90% profile values respectively as defined above. The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. The left and right field penumbra of the profile is given.

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to an inaccurate penumbra value being reported.

## Inflection Point Penumbra Parameters

The inflection point penumbra parameters are calculated relative to the dose value of the inflection point.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are available:

### Inflection Point Penumbra 80%-20%

The traditional parameters for defining the penumbra do not work for FFF beams. The 80%-20% penumbra based on the inflection point is calculated as the distance between the *80% left inflection point distance* to the *20% left inflection point distance* for the left penumbra and similarly for the right penumbra.



The profile must be centered on the origin as defined by the array or imaging modality for best results.

## Penumbra Slope

The traditional parameters for defining the penumbra do not work for FFF beams. The slope of the Hill function at the inflection point provides a good indication of the width of the penumbra. The *Inflection point left slope* and *Inflection point right slope* are given.



The profile must be centered on the origin as defined by the array or imaging modality for best results.

## Symmetry Parameters

Currently all symmetry parameters are calculated relative to the 50% edge or to the in field area as defined by the 50% edge.


## Area Symmetry

The area symmetry is the normalised difference between the right and left areas under the profile up to the *field edges*:

$$100 \cdot \frac{RA - LA}{RA + LA}$$

Where *RA* and *LA* are the area under the profile from the profile or detector centre to the right and left *Field Edge* respectively.



The area symmetry is affected by the field centre. If the field is slightly offset you can use the “Centre field” tool  to correct any offset.



Some programs use a multiplicative factor of 200.


## Maximum Dose Ratio

The maximum dose ratio symmetry (also called the point difference quotient) is the maximum ratio between the left and right profile values at the same distance from the profile or detector centre:

$$100 \cdot \text{maximum} \left[ \frac{P(dL)}{P(dR)}, \frac{P(dR)}{P(dL)} \right]$$

for  $dR = -dL$  from the origin as defined by the array or imaging modality to the end of the in field area.



The ratio symmetry is affected by the field centre. If the field is slightly offset you can use the “Centre field” tool  to correct any offset.

## Percent Symmetry

Percent symmetry is the ratio symmetry expressed as a percentage:

$$(RS - 1) * 100$$

where *RS* is the *Maximum Dose Ratio*.


## Point Difference

The point difference (also called maximum variation) is the maximum absolute difference between the left and right profile values at the same distance from the profile or detector centre:

$$100 \cdot \frac{P(dL) - P(dR)}{cax}$$

for  $dR = -dL$  from 0 to 80% of the field size



The point difference symmetry may be affected by the field centre. If the field is slightly offset you can use the “Centre field” tool  to correct any offset.

## Flatness Parameters

Flatness parameters are defined over the in field area.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are defined in the existing protocols:

### Dose Difference

The flatness dose difference parameter is the normalised difference between the maximum and minimum profile values taken over 80% of the field size:

$$100 \cdot \frac{max - min}{max + min}$$

Where *max* and *min* are the profile maximum and minimum respectively.



The *max* and *min* may be affected by the set window levels.

### Dose Ratio

Dose ratio flatness is the ratio of the maximum profile value to the minimum profile value expressed as a percent taken over the in field area:

$$100 \cdot \frac{max}{min}$$

Where *max* and *min* are the profile maximum and minimum respectively.



The *max* and *min* may be affected by the set window levels.

### Mean Value

The mean value dose parameter is the normalised average between the maximum and minimum profile values taken over 80% of the field size:

$$\frac{100}{cax} \cdot \frac{max + min}{2}$$

Where *max* and *min* are the profile maximum and minimum respectively and *cax* is the profile centre value.



The *max*, *min* and *cax* may be affected by the set window levels.

### Flatness L90/L50

The 90%/50% isodose ratio is defined as the ratio of the distance to the first point on the profile at 90% of the CAX value to the first point that reaches 50% of the CAX value or field edge. :

$$L90/L50 = maximum \left[ \frac{dl(P_{90\%})}{dl(P_{50\%})}, \frac{dr(P_{90\%})}{dr(P_{50\%})} \right]$$

$$P_{50\%} = CAX * 0.5$$



$$P_{90\%} = CAX * 0.9$$

Where  $dl(P_{50\%})$  and  $dl(P_{90\%})$  are the linear distances to the 50% and 90% profile values respectively on the left side of the profile as defined above, and  $dr(P_{50\%})$  and  $dr(P_{90\%})$  are the linear distances to the 50% and 90% profile values respectively on the right side of the profile as defined above. The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. The maximum ratio between the left and right sides is taken.

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to an inaccurate penumbra value being reported.

## Flatness L80/L50

The 80%/50% isodose ratio is defined as the ratio of the distance to the first point on the profile at 80% of the CAX value to the first point that reaches 50% of the CAX value or field edge. :

$$L80/L50 = maximum \left[ \frac{dl(P_{80\%})}{dl(P_{50\%})}, \frac{dr(P_{80\%})}{dr(P_{50\%})} \right]$$

$$P_{50\%} = CAX * 0.5$$

$$P_{80\%} = CAX * 0.8$$

Where  $dl(P_{50\%})$  and  $dl(P_{80\%})$  are the linear distances to the 50% and 80% profile values respectively on the left side of the profile as defined above, and  $dr(P_{50\%})$  and  $dr(P_{80\%})$  are the linear distances to the 50% and 80% profile values respectively on the right side of the profile as defined above. The algorithm searches symmetrically out from the profile centre on both sides of the profile to these points. The maximum ratio between the left and right sides is taken.

Since version 0.5 BeamScheme no longer does automatic grounding of the profile, i.e. reducing the minimum to zero. This means that if you want the grounded field size you must *explicitly* normalise or window the image.



Profiles are no longer grounded automatically. This can lead to differences in results from previous version.



Noisy data can lead to an inaccurate penumbra value being reported.

## Deviation Parameters

Deviation parameters generally apply only to WFF beams and express the amount of deviation from the CAX value over the in field area.

Parameters are calculated by combining base parameters in an arithmetic expression. The following parameters are defined in existing protocols:

### MAX/CAX

MAX/CAX flatness is the ratio of the CAX profile value to the maximum profile value expressed as a percent taken over 80% of the field size:

$$100 \cdot \frac{max}{cax}$$

Where *max* and *cax* are the profile maximum and central axis values respectively.



The *max* and *cax* are affected by the set window levels.

## Maximum Variation

This is the maximum variation in the ration of the absorbed dose at any point in the flattened area to that of the central axis expressed as a percentage difference between the lowest and highest value of this ratio.

$$100 \cdot \frac{max - min}{cax}$$

Practically though this has been implemented as:

$$100 \cdot maximum \left[ \frac{max - cax}{cax}, \frac{min - cax}{cax} \right]$$

Where *max*, *min* and *cax* are the profile maximum, minimum and central axis values respectively.



The *max*, *min* and *cax* are affected by the set window levels.

## 9.3 FFF algorithms

For FFF beams the standard interpolated or dose level parameters do not give good results. FFF fields are sharply peaked and the FWHM (50% dose level) is not a good indicator of the collimator field size. For FFF beams the inflection point (the point where the slope of the penumbra is greatest) of the penumbra has been proposed to represent the field edge.

While the maximum slope of the penumbra can be calculated directly this can be quite inaccurate for sparse profiles and recommended method is to fit a sigmoid or Hill function to the penumbra data and determine the inflection point of the function.

This is described under *Inflection Point Parameters*

## 9.4 In field area

The in field area (IFA) refers to the flattened area of WFF beams. A number of varying definitions exist depending on the protocol. Practically these are very difficult to implement as they are usually only defined for certain profile angles and field sizes. An additional problem is that the IFA is defined on the nominal field width which may be widely different from the actual field width.

To avoid these problem BeamScheme simply takes the in field area as 80% of the FWHM. This corresponds to most definitions for standard fields such as 10x10 and 20x20, but may overestimate the in field area for small fields (< 5cm) and underestimate the in field area for large fields (> 30cm).

The field statistics, flatness and symmetry values are all calculated from the in field area. For any discrepancies between values calculated by BeamScheme and other software first examine the defined in field area.

## PROTOCOLS

From version 0.5 BeamScheme supports definable protocol and parameter sets. A library of standard protocols is supplied. These may be edited and renamed or new protocols defined.

### 10.1 Selecting a protocol

Use the drop down list in the Results window to select the desired protocol. The parameters and results are displayed in the grid below.

### 10.2 Defined protocols

BeamScheme has a library of standard protocol sets. Do not alter these protocols. Instead copy them and create your own sets if you wish to make changes. There is considerable variation in the implementation of the protocols by different vendors. BeamScheme attempts to provide the consensus implementation, but there is no guarantee that this will match a particular vendor.

The available predefined protocols are:

#### 10.2.1 Default

This protocol gives the parameters as calculated by previous versions of BeamScheme.

The defined parameters are

- *Field Edge*
- *Field Centre*
- *Field Size*
- *Penumbra 90%-10%*
- *Penumbra 80%-20%*
- *Penumbra 90%-50%*
- *Area Symmetry*
- *Maximum Dose Ratio*
- *Percent Symmetry*
- *Dose Difference*

- *Dose Ratio*
- *MAX/CAX*

#### Flattened field

BeamScheme uses 80% of the *field size* as the flattened area.

### 10.2.2 AFFSAPS-JORF

#### Photon Profile Parameters

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Mean Value*
- *MAX/CAX*

#### Electron Profile Parameters

The defined electron parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Flatness L90/L50*
- *Maximum Variation*

#### Flattened field

BeamScheme uses 80% of the *field size* as the flattened area.

### 10.2.3 DIN

#### Photon Profile Parameters

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Dose Ratio*

**Electron Profile Parameters**

The electron profile parameters are the same as the photon profile parameters.

**Flattened field**

BeamScheme uses 80% of the *field size* as the flattened area.

## 10.2.4 IEC 60976

**Photon Profile Parameters**

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Mean Value*
- *Maximum Variation*

**Electron Profile Parameters**

The defined electron parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Flatness L90/L50*
- *Maximum Variation*

**Flattened field**

BeamScheme uses 80% of the *field size* as the flattened area.

## 10.2.5 Elekta

**Photon Profile Parameters**

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Mean Value*
- *Maximum Variation*

**Electron Profile Parameters**

The defined electron parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Maximum Dose Ratio*
- *Flatness L90/L50*
- *Maximum Variation*

#### **Flattened field**

BeamScheme uses 80% of the *field size* as the flattened area.

### **10.2.6 Siemens**

#### **Photon Profile Parameters**

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Area Symmetry*
- *Dose Difference*
- *MAX/CAX*

#### **Electron Profile Parameters**

The electron parameters are the same as the photon parameters

#### **Flattened field**

BeamScheme uses 80% of the *field size* as the flattened area.

### **10.2.7 Varian**

#### **Photon Profile Parameters**

The defined photon parameters are:

- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Point Difference*
- *Dose Difference*

#### **Electron Profile Parameters**

The electron parameters are the same as the photon parameters

#### **Flattened field**

BeamScheme uses 80% of the *field size* as the flattened area.

### 10.2.8 FFF

#### Photon Profile Parameters

FFF parameters only apply to photon beams:

- *Field Edge*
- *Field Centre*
- *Field Size*
- *Penumbra 80%-20%*
- *Penumbra Slope*
- *Area Symmetry*

#### Flattened field

BeamScheme uses 80% of the *field size* as the flattened area.





## ERRORS AND TROUBLE SHOOTING

BeamScheme has sophisticated error checking and detection. Error messages and warnings are displayed in the Status Bar. During the operation of the program you may encounter a variety of warning and error messages. There are three types of messages.

- 1) Informational: These give details of the programs working. These messages have a grey background. No action is required they are for information only. For example the startup message “BeamScheme initialised correctly”.
- 2) Actional: The program is requiring some kind of response from you. These messages have a yellow background. For example “Click Point A”
- 3) Error: The program has detected an error condition. These messages are displayed with a red background. Some errors will not allow you to continue until the error has been corrected. In this case the area or areas with the error will be highlighted in red.

Some common error messages and their causes are discussed below:

### 11.1 Unrecognised file or file read error

BeamScheme does not recognise this type of file. Please check that the file is in the list of *Supported File Formats*. Also check that the file does not contain compression or encryption.

### 11.2 File error, no data found

BeamScheme could not find readable image data in the file. The file may contain compression or be corrupt.

### 11.3 This software can not read compressed or 24-bit color files

The DICOM image is compressed or contains a colour image. Decompress the image using a PACS system such as CONQUEST or convert the colour image to monochrome using an editor such as the GIMP.

### 11.4 Unable to load DICOM header segment. Is this really a DICOM compliant file?

The file extension is “.dcm” but BeamScheme cannot find valid DICOM header tags. Try and open the file in another DICOM reader.

## 11.5 This file does not have enough data for the image size

The DICOM file size is less than the image dimensions multiplied by the pixel depth. Check that the file is not corrupted or compressed.

## 11.6 File error, corrupt file

A fatal error occurred while trying to read the file. Try and open it with another reader. If this works and you feel BeamScheme should be able to read this file please contact *Technical Assistance*.

## 11.7 Could not save protocol

An error occurred while trying to save the protocol. Check that you have rights to write to the directory and that there is enough disk space.

## 11.8 No protocol definition files found. Please create a file

You must have at least one protocol definition file available either in the 'Protocol' subdirectory of the executable (program) directory, or in the user data directory. On installation BeamScheme will copy a number of default protocols into the 'BeamScheme' subdirectory of the user configuration directory. Any new protocol is automatically saved in the 'BeamScheme' subdirectory of the user configuration directory.

Check that you have not accidentally deleted the protocols. If you have you will have to reinstall BeamScheme.

## 11.9 Could not read RAW text file

BeamScheme could not read values from a raw text image file. Values should be text and separated by spaces.

## 11.10 Could not evaluate expression

The protocol expression you have just entered contains a mistake. Expressions should not contain spaces or any special characters other than the standard pascal arithmetic operators.

The parser also may not have recognised the name of a function or variable you have used. Please check that the functions are valid pascal functions and that the variable names correspond to the defined variables.

## 11.11 Unable to fit curve! Initial parameter equal to zero

The non-linear Hill regression failed to fit the data. Check that the penumbra is well formed and that it does not cross the origin.

Other errors may from time to time be generated by the operating system or math subsystem. These will generally be in the form of a dialogue box requiring acknowledgement. Please refer to the appropriate documentation and/or contact technical assistance.

## TECHNICAL ASSISTANCE

There are, of course, no bugs in the software. However, it may happen that the software does not behave in a manner you expect or generates results that you think are wrong. If so, we would like to know about it. Likewise, if you feel BeamScheme should be doing something that it isn't, please tell us.

The authors can be contacted on [alanphys@yenzakahle.co.za](mailto:alanphys@yenzakahle.co.za).



## INDICES AND TABLES

- `genindex`
- `search`



## A

Add  
     Parameter, 15  
 Algorithms, 20  
     FFF, 38  
     In field area, 38  
 Angle  
     Profile, 16

## B

BeamScheme  
     Running, 9  
     Using, 11

## C

Centre field  
     Image, 18  
 Conventions, 3

## D

Delete  
     Parameter, 15  
 Deviation  
     MAX/CAX, 37  
     Maximum Variation, 38  
 DICOM  
     File format, 19  
 Disclaimer, 1

## E

Edit  
     Protocol, 14  
 Editing  
     Exit, 15  
     Protocol, 14  
 Errors, 43  
 Exit  
     Editing, 15  
     File, 13  
 Export  
     Menu, 13  
     X Profile, 13, 14

Y Profile, 14

## F

FFF  
     Algorithms, 38  
 Field  
     CAX Dose, 22  
     Centre, 31  
     Centre inflection point, 32  
     Edge, 31  
     Edge Inflection, 32  
     Maximum dose, 22  
     Maximum position, 23  
     Minimum dose, 23  
     Number of points, 23  
     Size, 31  
     Size inflection point, 32  
     Statistics, 22  
     Sum, 23  
     Sum squares, 23  
 File  
     Exit, 13  
     Menu, 13  
     Open, 13  
 File format  
     DICOM, 19  
     HIS, 20  
     IBA, 19  
     Image, 19  
     iPlan, 20  
     MapCheck, 19  
     PTW, 19  
     XiO, 20  
 File formats, 18  
 Flatness  
     Dose difference, 36  
     Dose ratio, 36  
     L80/L50, 37  
     L90/L50, 36  
     Mean value, 36

## H

Help

- Menu, 16
- HIS
  - File format, 20

## I

- IBA
  - File format, 19
- Image
  - Centre field, 18
  - File format, 19
  - Invert, 17
  - Normalise to CAX, 17
  - Normalise to MAX, 17
  - Pane, 17
- In field area
  - Algorithms, 38
- Installation, 7
- Installing
  - Linux, 9
  - Windows, 9
- Introduction, 1
- Invert
  - Image, 17
- iPlan
  - File format, 20

## L

- Licence, 1
- Linux
  - Installing, 9
  - Running, 11
  - Uninstalling, 9

## M

- MapCheck
  - File format, 19
- Menu, 13
  - Export, 13
  - File, 13
  - Help, 16
  - Print, 15
  - Window, 16

## N

- Normalise to CAX
  - Image, 17
- Normalise to MAX
  - Image, 17

## O

- Offset
  - Profile, 16
- Open
  - File, 13

## P

- Parameter
  - Add, 15
  - Delete, 15
- Parameters, 22
  - ADiff, 30
  - ALeft, 29
  - ARight, 29
  - Base, 22
  - calculated, 30
  - deviation, 37
  - Dose point, 28
  - field, 30
  - field inflection, 31
  - field interpolated, 31
  - Flatness, 35
  - IL20, 27
  - IL50, 27
  - IL80, 26
  - Inflection point, 25
  - Interpolated, 23
  - IR20, 27
  - IR50, 27
  - IR80, 26
  - L10, 25
  - L20, 25
  - L80, 24
  - L90, 24
  - LD20, 28
  - LD50, 28
  - LD60, 29
  - LD80, 29
  - Left edge, 24
  - LInf, 26
  - LSlope, 27
  - Penumbra, 32
  - Penumbra inflection point, 34
  - Penumbra interpolated, 32
  - R10, 25
  - R20, 25
  - R80, 24
  - R90, 24
  - RD20, 28
  - RD50, 28
  - RD60, 29
  - RD80, 29
  - RDiff, 30
  - Right edge, 24
  - RInf, 26
  - RSlope, 28
  - Symmetry, 29
  - symmetry, 34
- Penumbra
  - 80-20%, 33



- 80-20% inflection point, [34](#)
- 90-10%, [33](#)
- 90-50%, [33](#)
- slope, [34](#)
- Print
  - Menu, [15](#)
- Print PDF, [15](#)
- Profile
  - Angle, [16](#)
  - Offset, [16](#)
  - Pane, [16](#)
  - Width, [16](#)
- Profiles
  - diagonal, [21](#)
  - skew, [21](#)
- Protocol
  - AFFSAPS-JORF, [40](#)
  - Default, [39](#)
  - DIN, [40](#)
  - Edit, [14](#)
  - Editing, [14](#)
  - Elekta, [41](#)
  - FFF, [42](#)
  - IEC-60976, [41](#)
  - Save, [15](#)
  - Selecting, [39](#)
  - Siemens, [42](#)
  - Varian, [42](#)
- Protocols, [38](#)
  - defined, [39](#)
- PTW
  - File format, [19](#)

## R

- Requirements
  - System, [5](#)
- Results, [18](#)
- Running
  - BeamScheme, [9](#)
  - Linux, [11](#)
  - Windows, [11](#)

## S

- Save
  - Protocol, [15](#)
- Selecting
  - Protocol, [39](#)
- Symmetry
  - Area, [34](#)
  - Maximum dose ratio, [35](#)
  - Percent, [35](#)
  - Point difference, [35](#)
- System
  - Requirements, [5](#)

## T

- Technical Assistance, [46](#)
- Trouble shooting, [43](#)

## U

- Uninstalling, [9](#)
  - Linux, [9](#)
  - Windows, [9](#)
- Using
  - BeamScheme, [11](#)

## W

- Width
  - Profile, [16](#)
- Window
  - Menu, [16](#)
- Windows
  - Installing, [9](#)
  - Running, [11](#)
  - Uninstalling, [9](#)

## X

- X Profile
  - Export, [13](#), [14](#)
- XiO
  - File format, [20](#)

## Y

- Y Profile
  - Export, [14](#)